


Skin Grafts - From Macro to Minor

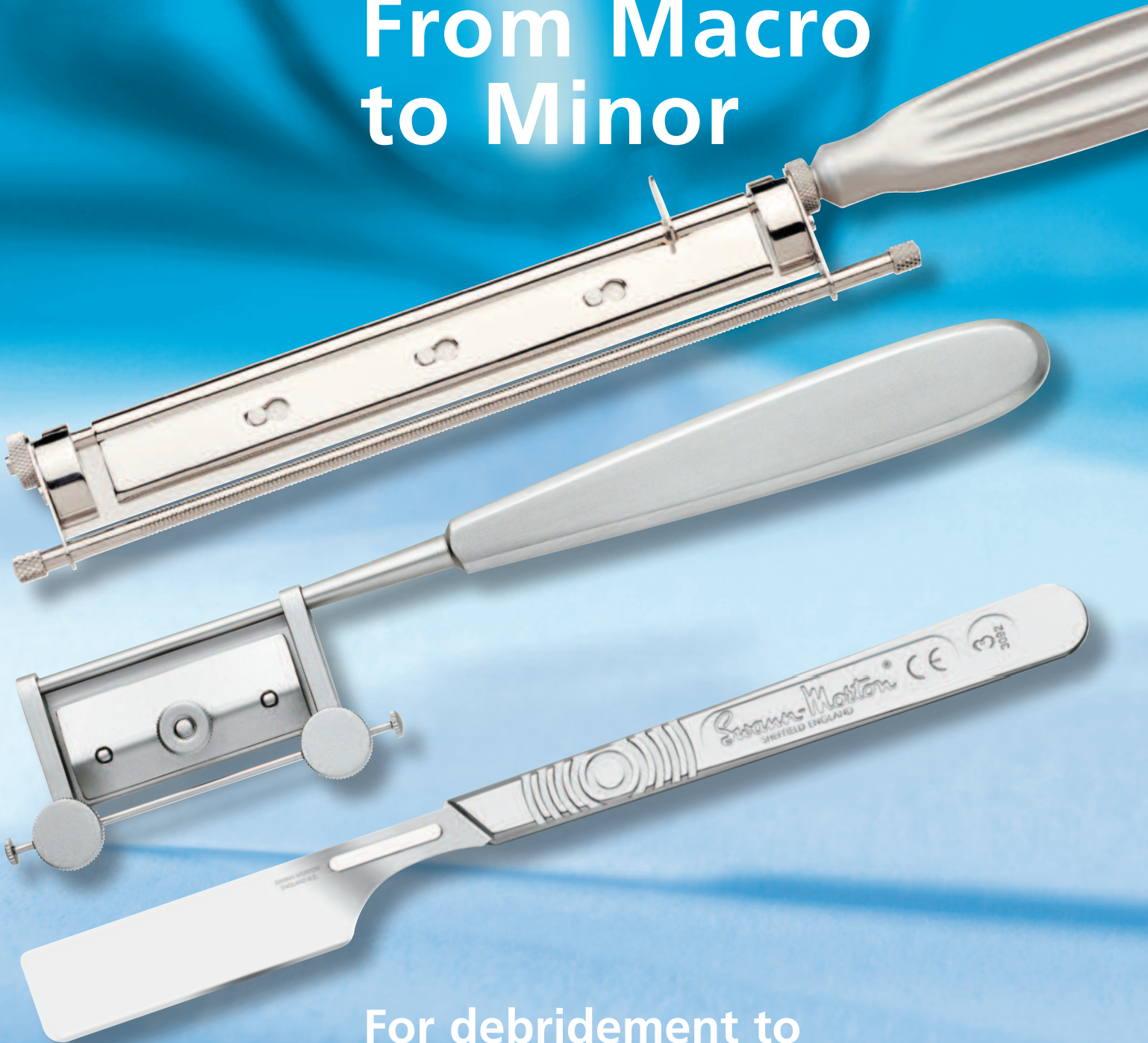

For debridement to
harvesting of skin

Swann-Morton®

The Skin Graft range is used in the surgical removal of a patch of healthy skin from one area of the patient's body which is then transplanted to another area where there has been extensive tissue loss. Skin grafts may be recommended for burns, specific surgeries that may require skin grafts for healing to occur and for areas of prior infection with extensive skin loss. They may also be used for cosmetic reasons during reconstructive surgeries.

With three models to choose from, namely the Braithwaite, Cobbett and Watson, the Stainless Steel Knife Handles are used in conjunction with the sterile 158mm Skin Graft Blades for harvesting "Macro" or "Major" Skin Grafts. The knives have an in-built adjustable guard which can limit the maximal thickness of the graft harvested depending on the properties of the skin at the donor site. This can vary in thickness due to multiple factors including the age of the patient.

The exact thickness of the graft can also depend upon:

- 1) The contact angle between the blade and skin whilst harvesting
- 2) The pressure applied by the operator onto the knife whilst harvesting

These are important aspects of skin grafting technique that need to be mastered to enable reproducibly consistent harvesting of high quality split thickness skin grafts.



The smaller Silvers stainless steel "miniature" Skin Graft Knife is available with either a distinctive blue autoclaveable handle or a solid stainless steel handle and was originally developed by H.L. Silver in Toronto, Canada.

The "new" custom designed surgical stainless steel Silver's blade has a 1½ inch (4 cm) postage stamp size cutting edge. The steel thickness now mirrors that of other skin graft blades such as the 9940 - 158 mm providing more stability for the user and preventing the occurrence of skin tear or roll up during harvesting. They are ideal for harvesting small intricate grafts from areas which would not be accessible with the larger standard knives such as the ear, dorsa of the hands or scalp. Although used mainly within Plastic, Reconstructive and Oral surgery the Silvers can also be used for the debridement of wounds.

The Dermatome blade range is used in the surgical removal of damaged tissue and harvesting of "borrowed" tissue from a donor site. The blades are fitted into a handheld instrument which resembles an electric razor.

The blade oscillates and the movement back and forth will "evenly" skim off the surface layers of skin. Electric Dermatomes are preferable for cutting thinner and longer strips of skin with a more homogeneous thickness.

Available in three sizes, a 25mm single segment (9943), a 50mm two segmental (9944) and a 75mm three segmental (9945), the blades are individually sealed in foil packets and boxed in units of 10 pieces. The blades are compatible with most leading brands of Dermatome knives.

The sterile SG3 minor skin graft blade can be conveniently fitted onto a standard No 3 surgical handle and was designed to assist with the debridement of burn related wounds at the site of major accidents and disasters.

It has been used on the hands and feet of small children given its practical size advantage compared to a standard skin grafting knife and also within Histology departments for the fine sectioning of the pathology blocks. It has a standard No 3 fitment so is compatible with handles 3 Graduated, 3L, 5B, 7, 9, B3 and B3L.

The latest addition to the burns and reconstructive surgery range is the SGD sterile single use Disposable scalpel.

This double edged minor skin graft blade, attached to a plastic handle for convenience, was originally introduced to assist in reducing the occurrence of post procedure infections amongst infants having undergone circumcision and has since become an essential component of procedure kits available at local pharmacy level. The scalpels are provided individually sealed in a tray and boxed in units of 10.

Following the popularity of the SGD sterile single use Disposable Scalpels a demand has arisen for the supply of the double edged stainless steel SGD blade to be used in conjunction with the standard Stainless Steel No.3 fitment surgical handle.

The SGD blades are provided sterile, individually sealed in foil packs and boxed in units of 20.

More information on the Swann-Morton Skin Graft Range can be found in the product section of our website at www.swann-morton.com.

Skin Graft Blades - Sterile

Blades in Stainless Steel, Sterile

Product	Product No	Unit
Minor SG3 to fit No 3 Handle		
SG3	2201	20 blades individually packed
SGD Blades - Sterile		
SGD	2203	20 blades individually packed
SGD Disposable Scalpel - Sterile		
SGD scalpel	0527	10 scalpels individually packed

Skin Graft Blades - Sterile

Blades in Stainless Steel, Sterile

Product	Product No	Unit
Silvers Skin Graft Blade for use with Silvers Handle		
Silvers	9942	10 blades individually packed
158mm Blade to fit Braithwaite, Cobbett, Watson		
158mm	9940	10 blades individually packed
Dermatome Blades - Sterile		
1 Segment	9943	10 blades individually packed
2 Segmental	9944	10 blades individually packed
3 Segmental	9945	10 blades individually packed

Skin Graft Handles - Non-Sterile


Stainless Steel, Non-Sterile, for use with Blade 9940

Product	Product No	Unit
Braithwaite Skin Graft Knife - Non-Sterile		
Skin Graft Knife	9901	1 handle individually boxed
Cobbett Skin Graft Knife - Non-Sterile		
Skin Graft Knife	9902	1 handle individually boxed
Watson Skin Graft Knife - Non-Sterile		
Skin Graft knife	9903	1 handle individually boxed

Skin Graft Sets

Stainless Steel Handle, 2 boxes of 10, sterile 158mm 9940 Blades

Product	Product No	Unit
Braithwaite Skin Graft Knife Set		
Skin Graft Knife Set	9904	1 Set
Cobbett Skin Graft Knife Set		
Skin Graft Knife Set	9905	1 Set
Watson Skin Graft Knife Set		
Skin Graft Knife Set	9906	1 Set

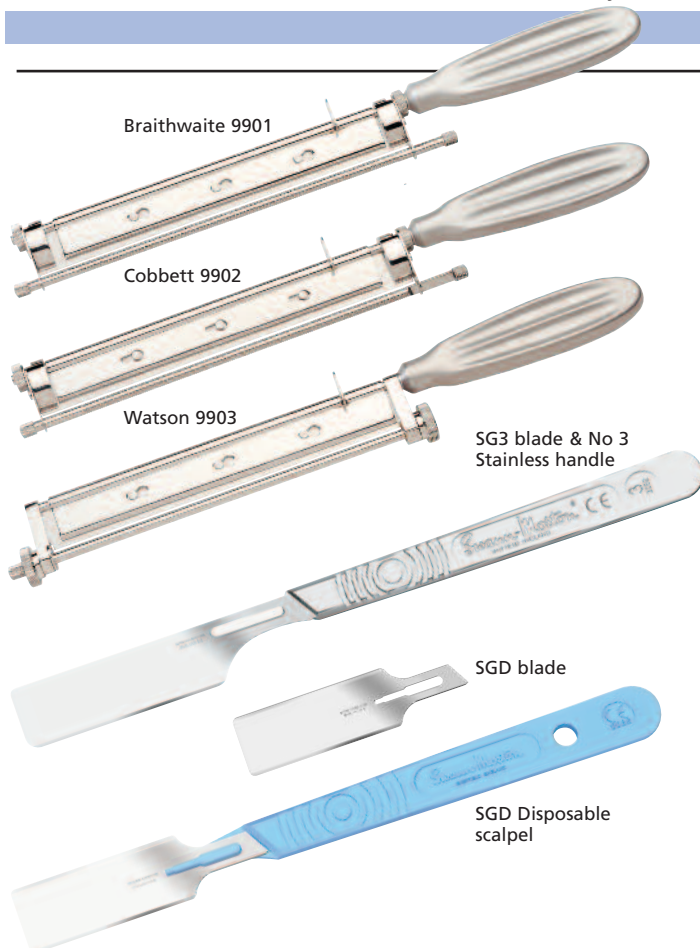

Skin Graft Handles - Non-Sterile

Stainless Steel Non-Sterile, for use with Blade 9942

Product	Product No	Unit
Silver Skin Graft Knife - Non-Sterile		
Silvers Skin Graft Knife	9911 (Blue)	1 handle individually boxed
Silvers Skin Graft Knife	9912 (Silver)	1 handle individually boxed



Skin Grafts - From Macro to Minor

A comprehensive range of Skin Graft products designed specifically for debridement to harvesting of skin.

The Skin Graft range is used in the surgical removal of a patch of healthy skin from one area of the patient's body which is then transplanted to another area where there has been extensive tissue loss. Skin grafts may be recommended for burns, specific surgeries that may require skin grafts for healing to occur and for areas of prior infection with extensive skin loss. They may also be used for cosmetic reasons during reconstructive surgeries.



Scan this code for more product information

Swann-Morton and the 'ring pattern' logo are registered trademarks of Swann-Morton Limited and related companies. The products shown are manufactured for use by trained health care professionals and can cause injury if used or handled incorrectly. The information and illustrations used in connection with the products shown on this leaflet are for general guidance only and users must seek proper training and instruction in their use before using or handling such products. Swann-Morton Limited exclude all liability for injury to the extent permitted by law. Products images for illustration purposes only and are not to scale.

Swann-Morton[®]
Sheffield - England

Swann-Morton Ltd. Penn Works, Owlerton Green, Sheffield S6 2BJ
Telephone: 0114 2344231, Sales 0114 2344223, Fax: 0114 2314966
uksales@swann-morton.com, exportsales@swann-morton.com
orders@swann-morton.com, www.swann-morton.com