

The redefined choice



The Swann-Morton Dermaplaning range - consistent quality, reliability and flexibility

Dermaplaning is a non-invasive aesthetics treatment which uses a surgical blade, scalpel or razor to carefully scrape the surface of the skin, removing unwanted hair, 'peach fuzz', also referred to as vellus hair from the epidermis in areas such as the upper lip, eyebrow, forehead, cheeks, and chin.

Dermaplaning can also :

- Remove dead skin cells from the face and neck.
- Treat the face and neck with light exfoliation.
- Reduce fine lines and wrinkles due to ageing, stress, smoking, and other physical or environmental factors.
- Diminish acne scars and visible pores, blotching, or other skin impurities.
- Support professional-looking brow shaping.

Swann-Morton®
SHEFFIELD ENGLAND



**UK
CA**
All medical devices comply with the Medical Device Directive 93/42/EEC as referred to in Part II of the Medical Devices Regulations 2002, Annex II excluding section 4 (as modified by Part 2 of Schedule 2A to the Medical Devices Regulations 2002)



Scan this code for more product information



ONE OF THE WORLD'S LEADING SURGICAL BLADE MANUFACTURERS

For Dermaplaning, sterile surgical blades and scalpels would traditionally be used by trained Aestheticians (Estheticians) in a controlled professional environment such as a Salon, Spa or Surgery. Quality surgical blades are often the number one choice of these trained aestheticians and cosmetically trained beauticians due to the consistent quality, reliability and flexibility of the blade and scalpel handle combination.

The most popular blade shapes for Dermaplaning are the Nos.10R and 14 which are also available as a complete scalpel with a blue handle already attached and even in a safety scalpel format with a high visibility red handle and retracting blade.



The 10R is often referred to as the 'butter blade'. It's short length and rounded tip making it ideal to treat small uneven areas of the face and in avoiding accidental cuts of the skin.


The 14 is a longer blade and is used to treat the entire surface of the face. For maximum effect, you can use the 14 blades in tandem with the 10R,

ensuring you cover both large and smaller areas of the face providing an optimum finish.

Dermaplaning can be used as part of or in preparation for other aesthetics treatments such as chemical peel, microdermabrasion or Fluorescent Pulsed Light procedures.

Product	Product No	Unit
Dermaplaning Blades - Sterile		
Blades in Stainless Steel - Sterile, individually V.P.I. protected in Metal Foil Packets		
Fits No 3 handle		
No10R Stainless Sterile Blades	0390	Box of 100
No14 Stainless Sterile Blades	0319	Box of 100
Ancillary Items - Non Sterile		
Swann-Morton Blade Remover units	5525	Box of 1
No 3 Graduated Stainless Surgical Handles	0933	Box of 1

Product	Product No	Unit
Dermaplaning Scalpels - Sterile		
Blades in Stainless Steel - Sterile, individually V.P.I. protected in Metal Foil Packets		
No 3 handle		
No10R Sterile Retractable Scalpels	3990	Box of 25
No10R Sterile Retractable Scalpels	4990	Box of 10
No10R Sterile Disposable Scalpels with Guard	6690	Box of 10
No14 Sterile Disposable Scalpels	0519	Box of 10
No14 Sterile Retractable Scalpels	3919	Box of 25
No14 Sterile Retractable Scalpels	4919	Box of 10

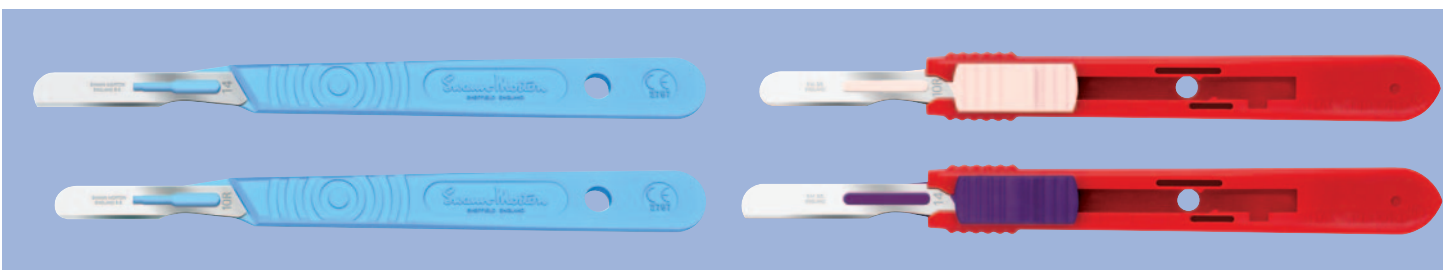

Swann-Morton and the 'ring pattern' logo are registered trademarks of Swann-Morton Limited and related companies. The products shown are manufactured for use by trained health care professionals and can cause injury if used or handled incorrectly. The information and illustrations used in connection with the products shown on this leaflet are for general guidance only and users must seek proper training and instruction in their use before using or handling such products. Swann-Morton Limited exclude all liability for injury to the extent permitted by law.